Tecnovita
by **BH**

Instrucciones de uso y garantía
Operating instructions and guarantee
Mode d'emploi et garantie
Instruções de utilização e garantia
Manuale d'uso e garanzia
Bedienungsanleitung und Garantie

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902 170 258

Thank you for purchasing our e-SCULPTER YR50 body sculptor.
All of our products undergo rigorous testing to ensure maximum quality.
To use the unit correctly, first read through the whole manual before using it.

INTRODUCTION


EMS, electric muscular stimulation is a widely recognized method that has been used in the sport field for a long time. The e-SCULPTER is both easy and comfortable to use. This product has been designed specifically to help tone abdominal, legs and arms muscles. It allows you to enjoy the exercise and increase muscle capacity by developing strength and improving tone effortlessly. Thanks to its design, the e-SCULPTER tones muscles anytime, anywhere.

GENERAL SAFETY INSTRUCTIONS

1. Before using the appliance make sure it is in good working order. If you detect any faults get in touch with your local dealer.
2. The unit must not be used by children or the disabled without suitable adult supervision.
3. Pregnant women must not use it.
4. Do not use if you suffer from heart disease or other type of disease.
5. Do not use the appliance if you have cancer.
6. Do not use if you are prone to attacks, for example, epilepsy.
7. Do not use the appliance after surgical operations in which increased muscle tension may affect the recovery process
8. Do not use directly on open wounds or varicose veins.
9. Do not use if you have any implanted electrical device, such as a pacemaker.
10. Do not use if you have any metal implants.
11. Do not use near the genitals.
12. Do not use on broken or inflamed skin.
13. Do not use the appliance near your heart.
14. If you suffer from any serious illness, especially if you suspect that you may have high blood pressure or problems with blood clotting, consult your doctor before using the appliance.
15. Consult your doctor before using this product if you have any type of skin disease.
16. If you suffer from chronic pain of unknown origin, regardless of where on your body, consult your doctor before using the appliance.
17. People undergoing medical treatment should consult a doctor before using this appliance.
18. Diabetics should consult their doctor before using this appliance.
19. People with constant skin irritation resulting from prolonged stimulation by the same electrode, should consult a doctor.
20. Consult your doctor if you have any problems resulting from the stimulation treatment.
21. Stop using the appliance if you experience any discomfort during use.
22. Do not use after consuming alcohol.
23. If the appliance malfunctions or if you feel discomfort or pain, discontinue use immediately.
24. Do not use while sleeping, driving or while using some other device.
25. Keep the appliance away from heat sources and do not use close (<1 m) to medium or microwave wave devices (e.g. mobile phones) because they can cause power surges.
26. Do not use when there are high frequency surgical devices connected simultaneously.
27. Do not allow the electrodes to come into contact with metal objects, e.g. belts or collars, during operation.
28. Do not use near flammable, gaseous or explosive substances.
29. Switch the appliance off before removing it to prevent unnecessary irritation.

30. Never use the appliance without first applying plenty of electro-stimulation gel. If you do, then the stimulating effect will not be transmitted correctly and it can also damage the appliance.
31. Do not apply enriching creams beforehand, as these cause the electrodes to wear out faster or may cause unwanted current surges.
32. Do not use the unit in damp/humid surroundings and prevent it from getting wet. Do not touch the device with wet hands while it is operating. Under no circumstances must this product be used in the bath or shower.
33. Do not expose the appliance to direct sunlight or high temperatures.
34. Do not allow the appliance to come into contact with water or any other liquid. Do not immerse the appliance in liquid to clean it.
35. If you notice that the appliance is not operating correctly, do not attempt to repair, dismantle or modify it. Unauthorised repairs may cause injury, in addition to breaking the terms of the warranty. Contact your technical assistance service.

CONTENT

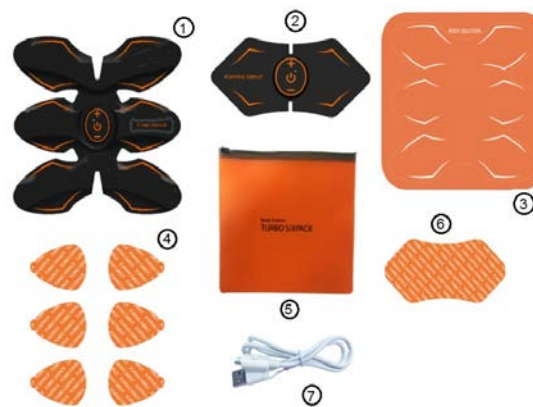1. SIXPACK SCULPTOR
2. BODY SCULPTOR
3. E-SCULPTOR HOLDER
4. SIXPACK GEL PADS
5. PVC BAG
6. BODY SCULPTOR GEL PADS
7. USB CABLE

DESCRIPTION

- Power: DC 5V;
- Output voltage: 55V a 500 Ω;
- Frequency: 0-100 Hz;
- Measurements: 18.5 x 17 cm and 14 x 6.5 cm;
- 5 pulse modes;
- 15 intensity levels;
- Countdown timer of 20 min.

Environmental conditions for operating

Temperature range: -10 / 40 °C

Humidity: 30-85% RH

OPERATING INSTRUCTIONS

1. Before use, it is recommended to charge the device with the enclosed USB cable.
2. Peel off the orange plastic film from one side of the enclosed electrode gel pad and install the electrode pad onto the back side of the control unit as shown in the picture.
3. Place the device onto the treatment areas, press down firmly and ensure a full and firm contact with skin.
4. Press the On/Off button to turn on the unit, indicated by the LED light and buzz sound for 1 second.
5. When turned on, the device works at combination mode.
6. Press the “+” button to increase the pulse intensity, and press the “-” button to decrease it. Each time the intensity level is changed by pressing “+” or “-”, the light will flash once and also the sound will buzz once.
7. Press the “+” button for 3 seconds to select a desired pulse mode according to the following order: 1,2,3,4,5. Similarly, pressing the “-” button for 3 seconds the selection will be done downward. (5,4,3,2,1). PLEASE NOTE: The Mode 1 is a combination of the others 4 modes. Each time the mode is changes, light and buzz sound will activate twice.
8. The countdown timer is 20 min. When time is over, the unit will turn off automatically. The unit could be also turned off anytime by pressing the On/Off button, indicated by a triple light flashing and a triple buzz sound.

Note: The device will automatically power off when not in use on any treatment area, indicating a light and a sound on for 3 seconds and then flashing and buzzing 3 times.

RECOMMENDED PRACTICE

- Treat each part of the body for a maximum period of 20 minutes.
- Start from the lowest intensity and gradually adjust to a comfortable level on a scale level from 1 to 15 making sure of a constant comfortable filling.
- Be sure that the area to be treated is free of perspiration, dirt, creams and oils and that the skin is healthy and free from abrasions.
- In order to extend their lifespan, keep electrode gel pads in the storage holder immediately after use.
- The electrode gel pads are disposable and should be replaced when they lose their adhesiveness.

TROUBLESHOOTING

PROBLEM	CAUSE	SOLUTION
Stimulation is weak or absent	-No contact with skin -Skin is not completely clean from creams, oils, etc. -Low battery	-Press the device onto the skin and ensure a better adhesiveness -Clean the skin from creams, oils, etc. -Charge the battery
The device does not turn on	- Low battery	- Charge the battery
The gel pads are not sticky anymore	-Lost of adhesiveness	-Replace the pads

MAINTENANCE

Turn off the device before any maintenance action.

It's recommended to clean the device with a smooth wet cloth.

Never use cleaning products or excessive water.

Gel pads remain stucked onto the device between uses and plastic film can be used to protect them from getting dirty and/or dry.

Remove gel pads only for its replacement.

Warranty

This product is guaranteed for 2 years in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25th May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10th July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT nº B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product.

Date of purchase.-	Purchased in.-
Purchaser.-	Dealer stamp.-

Please keep the receipt as proof of purchase.

Should you have any queries please do not hesitate to phone our customer service line:

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BH China Co., Ltd.